# RENPHO

---

User Manual

## Eye Massager

Model: R-G001

## ATTENTION

---

### READ THOSE INSTRUCTION

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, semental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
3. The battery must be removed from the appliance before it is scrapped.
4. The appliance must be disconnected from the supply mains when removing the battery.
5. The battery is to be disposed of safely.
6. WARNING: For the purposes of recharging the battery. only use the detachable supply unit provided with this appliance.
7. The appliance is only to be used with the power supply unit provided with the appliance.

8. This appliance contains batteries that are only replaceable by skilled persons.

9. This appliance is not intended for use by persons (including children ) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

**Please read this instruction manual carefully before using. RENPHO assumes no responsibility for any consequences caused by misuse or failure to follow these instructions.**

1. DO NOT use if you have had an eye operation, or if you have any eye diseases such as glaucoma, cataract, and amotio retinae.

2. Disabled people, and the elderly should only use this product when properly supervised.

3. DO NOT use if you have any swelling in the eye area.

4. If you feel any pain around your eyes, please discontinue using the device and consult a physician.

5. **DO NOT** use in high temperatures or humid environments. Avoid direct exposure to the sun.
6. **DO NOT** use the device with wet hands.
7. Remove contact lenses, glasses, false eyelashes, hairpins, and hair accessories, etc before using.
8. Only use the charger (5V 1A, not included) to recharge the eye massager.
9. **DO NOT** clean the product with water; clean it with cotton soaked in alcohol.
10. Avoid contact with sharp objects.
11. **DO NOT** use this product while driving.
12. If any fault or malfunction should occur during use, take it off immediately. Never attempt to disassemble or repair the device yourself.
13. The device is for private use only and cannot be used as medical equipment. It is not intended for medical diagnosis, treatment, or commercial purposes.
14. Use it in the morning and evening for no more than 10 minutes each time.
15. Please keep this manual for future reference.

# PRODUCT STRUCTURE

---

Power Button

Long Press: Power ON/OFF  
Short Press: Switch Modes




Mode Display

Adjustable Headband

Massage Area

Charging Port

USB Cable
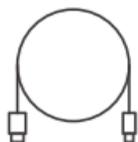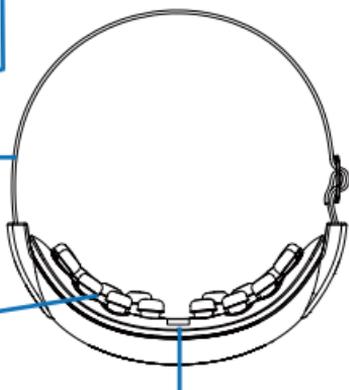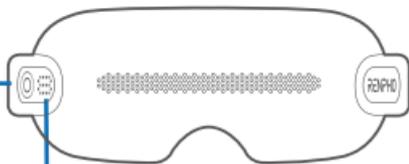
## SPECIFICATION

---

Product Name: Eye Massager

Model: R-G001

Input: 5V  1A

Rated Power: 5W

Default Auto-off Timer: 10 minutes

Battery Capacity: 1050mAh

Product Weight: 230g

Product Size: L196×W78×H82mm

Operation Frequency: 2402-2480MHz

Transmit Power: 8.81dBm MAX

## WHAT'S INCLUDED

---

1 x RENPHO Eye Massager

1 x USB Cable

1 x User Manual

## USING THE PRODUCT

---

### Charging

1. Use the 5V charger (not included) to charge your eye massager.
2. Insert USB cable into the socket of the mobile charger.
3. Softly insert the USB plug into the socket of the device,

- making sure the USB cable is inserted in the right direction.
4. While charging, the LED indicator will flash red. When done charging, the LED indicator will be purple.
  5. The charge duration will be about 2.5 hours. Disconnect the power adapter after charging is complete. Once fully charged, the massager can be used in 10-minute intervals up to 8 times.

- ※ **The massager and the adapter may become warm during normal charging process.**
- ※ **Use time may differ depending on the operation mode.**
- ※ **The device will not work while charging.**

### **Wear this device**

1. Loosen the adjustable headband.
  2. Place the device around your head and tighten the headband to the appropriate fit.
- ※ Be careful to ensure the device is not too tight around your head.
  - ※ Please make sure that you have removed your contact lenses, glasses, false eyelashes, etc. before using this device.
  - ※ Do not operate the device with wet hands.

## Bluetooth Connection

1. Turn on the device by long press the Power Button.
2. Enable your mobile phone Bluetooth setting and search for "R-G001", an audio voice will confirm when the device is connected to Bluetooth.
3. After successfully pairing your device, music can be pushed to the device for playing while massaging.

## Indicator Introduction

Status	Indicator Display
Normal Operation	Solid Purple
Low Battery	Flashing Red
Charging	Flashing Red
Fully Charged	Solid Purple
Over-voltage Charging (over 5V)	No light

## Operation

1. Long press the Power Button to turn on the device. The message will begin after the voice prompt is heard.
2. Short press the Power Button to switch between modes. An audio voice will confirm mode change and an LED light will display the current mode.

Mode 1: Tension Relief Mode

Mode 2: Retinal Protection Mode

Mode 3: Youthful Radiance Mode

Mode 4: Regeneration Mode

Mode 5: Rhythmic Mode (Vibration sync to the music beat)

3. Long press the Power Button to stop using the device.

4. A voice prompt will be heard when the device has been in operation for 10 minutes. An indicator light will go off and the device will automatically shut down.

### **Warm Tips**

※ When using this device, it's recommended to keep your eyes closed, lean back in the chair and completely relax your body.

※ Your vision will become blurry for a short period of time after use. Regular vision will return faster if eyes are kept closed for a short while.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Power-on failure	<ol style="list-style-type: none"><li>1. Check whether the battery level is too low</li><li>2. The time of pressing Power Button is too short</li></ol>	<ol style="list-style-type: none"><li>1. Please charge the battery timely</li><li>2. Long press the Power Button for 2-3 seconds</li></ol>
Low massage intensity	<ol style="list-style-type: none"><li>1. Battery level is low</li><li>2. The headband is too loose</li></ol>	<ol style="list-style-type: none"><li>1. Charge the device</li><li>2. Wear it tighter</li></ol>
Device turns off automatically	<ol style="list-style-type: none"><li>1. Battery is Low</li><li>2. 10 Minutes Automatic Shut-off</li></ol>	<ol style="list-style-type: none"><li>1. Charge the device</li><li>2. Turn on the massager again</li></ol>
The product can't be charged	<ol style="list-style-type: none"><li>1. The charge's voltage is too high</li><li>2. Adapter is broken</li></ol>	<ol style="list-style-type: none"><li>1. Use 5V 1A/2A charger to charge</li><li>2. Change another adapter</li></ol>
If you encounter other problems, turn off the device and contact us at <a href="mailto:support@renpho.com">support@renpho.com</a>		

## CLEANING AND MAINTENANCE

1. Please turn off the power when cleaning, do not use corrosive detergents such as benzene or diluents to clean the product.
  2. Do not expose the device to excessive sun or hot temperature.
- When the device has not been used for over 3 months, recharge it before use.

3. Do not place it next to a wall or on the ground or rooftop.
4. Keep the device dry, clean, and away from rats.
5. Do not place anything heavy on it. Always bring your product to a qualified professional for any service or repair.
6. Do not disassemble the product yourself.
7. Do not disregard the device with normal household waste at the end of its life; instead, follow local environmental regulations and guidance for disposal.

## WARRANTY (1-YEAR PERIOD)

The manufacturer's warranty is void under the following circumstances but not limited to:

- Improper or inadequate maintenance or modification.
- Accident, misuse, abuse, contamination, or other external causes.
- The use of unspecified adapter and accessories.
- Loss or damage in transit.
- Damages that occurs as a result of one's failure to follow the instructions.
- This warranty does not apply to expendable or consumable parts and does not extend to any product from which the serial number has been removed.

## CUSTOMER SERVICE

---

We have an award-winning warranty, exchange, and customer service team that guarantees hassle-free solutions to any issue you might have within the 24 business hours.

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 **Web: [www.renpho.com](http://www.renpho.com)**

## Notice:

This symbol on the device or its packaging indicates: Do not dispose of electrical appliances as unsorted municipal waste; use separate collection facilities. Contact your local authority for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. When replacing old appliances with new ones, the retailer is legally obliged to take back your old appliance for disposal at least free of charge.

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## Declaration of Conformity

We (EU Importer is responsible for this declaration)

C&EConnectionE-Commerce(DE)GmbH

(Company name)

ZumLinnegraben20,65933,FrankfurtamMain,Germany

(Company address)

declare under our sole responsibility that the product

Tradename: RENPHO

Equipment: Eye Massager

Model No. : R-G001

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to Radio Equipment Directive(2014/53/EU)& RoHS(2011/65/EU) product is responsible to affix **CE** marking, the following standards were applied:

*EN 62368-1:2014+A11:2017*

*EN 62479:2010*

*ETSI EN 301489-1 V2.2.3*

*ETSI EN 301489-17 V3.1.1*

*ETSI EN 300328 V2.2.2*

*IEC 62321 Series*

Full Name: Jacky

Position: Approbation Specialist

Signature:

Date: 2021-03-22

## Declaration of Conformity

We (EU Importer is responsible for this declaration)

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Datum: 2021-3-22

## Declaration of Conformity

We (UK Importer is responsible for this declaration)

Richers & Vision Ltd.

(Company name)

Unit 58B, Bordon Enterprise Park, Budds Lane, Bordon, Hampshire, GU35 0JE, U.K.

(Company address)

declare under our sole responsibility that the product

Tradename: RENPHO

Equipment: Eye Massager

Model No. : R-G001

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to Radio Equipment Directive(2014/53/EU)& RoHS(2011/65/EU) product is responsible to affix **CE** marking, the following standards were applied:

*EN 62368-1:2014+A11:2017*

*EN 62479:2010*

*ETSI EN 301489-1 V2.2.3*

*ETSI EN 301489-17 V3.1.1*


*ETSI EN 300328 V2.2.2*

*IEC 62321 Series*

Full Name: Zeke

Position: Approbation Specialist

Signature:

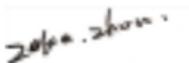Date: 2021-03-22

**Shenzhen Runfeng Digital Technology Co., LTD**

Chuanghui Building, Intersection of Banxuegang Boulevard and  
Yongxiang Road, Bantian Subdistrict, Longgang District, Shenzhen,  
China

**EU|REP| C&E Connection E-Commerce (DE) GmbH**

Zum Linnegraben 20, 65933, Frankfurt am Main, Germany  
Email: [info@ce-connection.de](mailto:info@ce-connection.de)

**UK|REP| Richers & Vision Ltd.**

Unit 58B, Bordon Enterprise Park, Budds Lane, Bordon, Hampshire,  
GU35 0JE, U.K.

Made in China  
**NA1**