RENPHO

User Manual

Eye Massager

Model: R-G001

ATTENTION

READ THOSE INSTRUCTION

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
3. The battery must be removed from the appliance before it is scrapped.
4. The appliance must be disconnected from the supply mains when removing the battery.
5. The battery is to be disposed of safely.
6. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
7. This appliance contains batteries that are only replaceable by skilled persons.

8. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

Please read this instruction manual carefully before using. RENPHO assumes no responsibility for any consequences caused by misuse or failure to follow these instructions.

1. DO NOT use if you have had an eye operation, or if you have any eye diseases such as glaucoma, cataract, and amotio retinae.
2. Disabled people, and the elderly should only use this product when properly supervised.
3. DO NOT use if you have any swelling in the eye area.
4. If you feel any pain around your eyes, please discontinue using the device and consult a physician.
5. DO NOT use in high temperatures or humid environments. Avoid direct exposure to the sun.
6. DO NOT use the device with wet hands.

7. Remove contact lenses, glasses, false eyelashes, hairpins, and hair accessories, etc before using.
8. Only use the charger (5V 1A, not included) to recharge the eye massager.
9. DO NOT clean the product with water; clean it with cotton soaked in alcohol.
10. Avoid contact with sharp objects.
11. DO NOT use this product while driving.
12. If any fault or malfunction should occur during use, take it off immediately. Never attempt to disassemble or repair the device yourself.
13. The device is for private use only and cannot be used as medical equipment. It is not intended for medical diagnosis, treatment, or commercial purposes.
14. Use it in the morning and evening for no more than 10 minutes each time.
15. Please keep this manual for future reference.

PRODUCT STRUCTURE

Power Button

Long press: Power ON/OFF

Short press once to switch modes.

Short press twice to switch between Silent Mode and Voice Mode.



Mode Display

Adjustable Headband

Massage Area

Charging Port

USB Cable

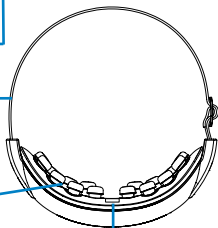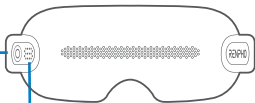SPECIFICATION

Product Name: Eye Massager

Model: R-G001

Input: 5V $\overline{=}$ 1A

Rated Power: 5W

Default Auto-off Timer: 10 minutes

Battery Capacity: 1050mAh

Product Weight: 230g

Product Size: L196×W78×H82mm

Operation Frequency: 2402-2480MHz

Transmit Power: 8.81dBm MAX

WHAT'S INCLUDED

1 x RENPHO Eye Massager

1 x USB Cable

1 x User Manual

USING THE PRODUCT

Charging

1. Use the 5V charger (not included) to charge your eye massager.
2. Insert USB cable into the socket of the mobile charger.
3. Softly insert the USB plug into the socket of the device, making sure

the USB cable is inserted in the right direction.

4. While charging, the LED indicator will flash red. When done charging, the LED indicator will be purple.

5. The charge duration will be about 2.5 hours. Disconnect the power adapter after charging is complete. Once fully charged, the massager can be used in 10-minute intervals up to 8 times.

The massager and the adapter may become warm during normal charging process.

Use time may differ depending on the operation mode.

The device will not work while charging.

Wear this device

1. Loosen the adjustable headband.

2. Place the device around your head and tighten the headband to the appropriate fit.

Be careful to ensure the device is not too tight around your head.

Please make sure that you have removed your contact lenses, glasses, false eyelashes, etc. before using this device.

Do not operate the device with wet hands.

Bluetooth Connection

1. Turn on the device by long press the Power Button.
2. Enable your mobile phone Bluetooth setting and search for "R-G001", an audio voice will confirm when the device is connected to Bluetooth.
3. After successfully pairing your device, music can be pushed to the device for playing while massaging.

Indicator Introduction

Status	Indicator Display
Normal Operation	Solid Purple
Low Battery	Flashing Red
Charging	Flashing Red
Fully Charged	Solid Purple
Over-voltage Charging (over 5V)	No light

Operation

1. Long press the Power Button to turn on the device. The massage will begin after the voice prompt is heard.
2. Short press the Power Button to switch between modes. An audio voice will confirm mode change and an LED light will display the current mode.

Mode 1: Tension Relief Mode

Mode 2: Retinal Protection Mode

Mode 3: Youthful Radiance Mode

Mode 4: Regeneration Mode

Mode 5: Rhythmic Mode (Vibration sync to the music beat)

3. Long press the Power Button to stop using the device.

4. A voice prompt will be heard when the device has been in operation for 10 minutes. An indicator light will go off and the device will automatically shut down.

Silent Mode & Voice Mode

Quick press the Power button twice to enter Silent Mode (without voice broadcast).

In the silent mode, quick press the Power button twice to switch back to the Voice Mode.

Warm Tips

When using this device, it's recommended to keep your eyes closed, lean back in the chair and completely relax your body.

Your vision will become blurry for a short period of time after use. Regular vision will return faster if eyes are kept closed for a short while.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Power-on failure	<ol style="list-style-type: none">1. Check whether the battery level is too low2. The time of pressing Power Button is too short	<ol style="list-style-type: none">1. Please charge the battery timely2. Long press the Power Button for 2-3 seconds
Low massage intensity	<ol style="list-style-type: none">1. Battery level is low2. The headband is too loose	<ol style="list-style-type: none">1. Charge the device2. Wear it tighter
Device turns off automatically	<ol style="list-style-type: none">1. Battery is Low2. 10 Minutes Automatic Shut-off	<ol style="list-style-type: none">1. Charge the device2. Turn on the massager again
The product can't be charged	<ol style="list-style-type: none">1. The charge's voltage is too high2. Adapter is broken	<ol style="list-style-type: none">1. Use 5V 1A/2A charger to charge2. Change another adapter
If you encounter other problems, turn off the device and contact us at support@renpho.com		

CLEANING AND MAINTENANCE

1. Please turn off the power when cleaning, do not use corrosive detergents such as benzene or diluents to clean the product.

2. Do not expose the device to excessive sun or hot temperature.

When the device has not been used for over 3 months, recharge it before use.

3. Do not place it next to a wall or on the ground or rooftop.
4. Keep the device dry, clean, and away from rats.
5. Do not place anything heavy on it. Always bring your product to a qualified professional for any service or repair.
6. Do not disassemble the product yourself.
7. Do not disregard the device with normal household waste at the end of its life; instead, follow local environmental regulations and guidance for disposal.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel:** +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. **DO NOT** dispose of any product parts as they may be required for inspection/repair.

Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.
EU REP C&E Connection E-Commerce (DE) GmbH
UK REP Richers & Vision Ltd.