
# RENPHO User Manual

## Leg Massager with Heat

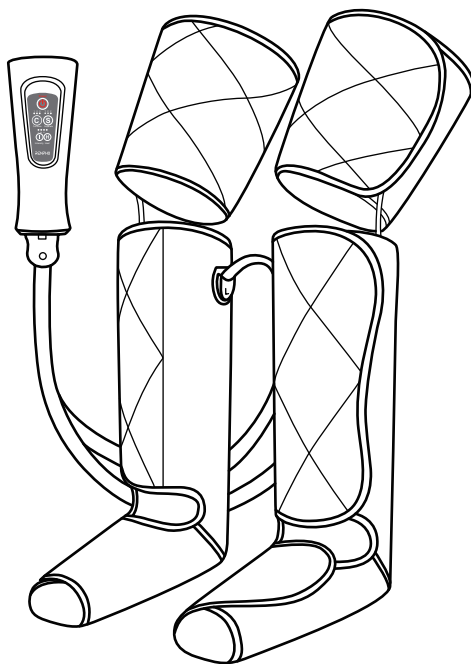Model: RP-ALM070H

# Contents

**English** 01~12

---

**Deutsch** 13~25

---

**Français** 26~38

---

**Italiano** 38~51

---

**Español** 52~64

---

Thank you for your purchase of RENPHO Leg Massager with Heat. This Leg Massager is designed to massage your feet, calves and thighs. It helps to relieve muscle pain and fatigue. You should use it with care and safety in mind. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

## Important Safety Instructions

When using the appliance, basic precautions should always be followed, including the following:

**READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.**

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or cognitive capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
2. Children should not play with the appliance.
3. Cleaning and user maintenance shall not be made by children.
4. This appliance should not be used for any medical purposes. It is designed for household use only.

**DANGER** - To reduce the risk of electric shock and fire:

1. Never use pins or other metallic fasteners with this appliance.
2. Ensure that the electrical voltage and frequency of the circuit correspond to those indicated on the adapter.
3. Do not leave the appliance unattended when it is working.
4. Always unplug the power adapter from the electrical outlet after each use.
5. Do not plug or unplug the power adapter or air hose with wet hands.
6. Keep the cord and appliance away from heated surfaces.
7. The appliance has a heating function. Persons insensitive/sensitive to heat should turn off the temperature while using.

**WARNING** -To reduce the risk of fire, electric shock, or injury to persons, please follow the instructions below:

1. If you fall into any of the following user groups, please consult a doctor before use:
  - Pregnant or lying-in women
  - Children
  - Elderly (over 65 years old)

- Persons with skin infections, sensitivity, redness, wounds, and/or burns on the skin
  - Persons using a pacemaker or other implantable medical device
  - Persons with heart disease, malignant tumors, or other serious health conditions
  - Persons receiving medical treatment or feel discomfort during massages
  - Persons suffering from osteoporosis, diabetes, or thrombosis
  - Persons with high or abnormal blood pressure
  - Persons suffering from a sensory impairment that reduces the feeling of pain or heat (e.g. metabolic disorder)- Persons suffering from fever over 100.4
2. Stop using if you feel abnormal or unwell after use and consult a doctor.
  3. **DO NOT** use spare parts that are not approved by RENPHO.
  4. Non-professional repair personnel should not dismantle or replace the appliance without manufacturer permission.
  5. **DO NOT** use the appliance if wrap surfaces or the power cord are damaged.
  6. **DO NOT** use or store the appliance in a place with humidity such as a bathroom or places with high temperature.
  7. **DO NOT** carry the appliance by the power cord or handheld controller.
  8. **DO NOT** use while sleeping, operating machinery, or after alcohol consumption.
  9. **DO NOT** stand or walk while using the appliance.
  10. **DO NOT** expose the handheld controller, the wraps, or the power adapter to water as this may result in malfunction or electric shock.
  11. **DO NOT** scratch, excessively bend, pull, or twist the wire of the power adapter.
  12. **DO NOT** place heavy objects on or insert any foreign objects into the slit (this may cause fire and electric shock).
  13. **DO NOT** use the sharp tools to scratch or puncture the wraps, as this may cause the appliance to an air leak.
  14. **DO NOT** use wraps directly on swollen, inflamed areas, sensitive skin, or skin eruptions.
  15. **DO NOT** use the appliance for more than 20 minutes at a time. Extensive use could lead to the appliance's excessive heating and shorter lifespan.
  16. **DO NOT** operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
  17. **DO NOT** use the appliance on pets or animals.
  18. Use heated surfaces carefully. May cause serious burns. The unattended use of heat by children or incapacitated persons may be dangerous.
  19. Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use it on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.

**SAVE THESE INSTRUCTIONS**


# About the RENPHO Leg Massager with Heat

## 1. What's in the box

1 x Pair of calf & foot wraps  
1 x Pair of thigh wraps  
1 x Handheld controller

1 x Power adapter  
1 x User manual  
1 x Storage bag


## 2. Product introduction

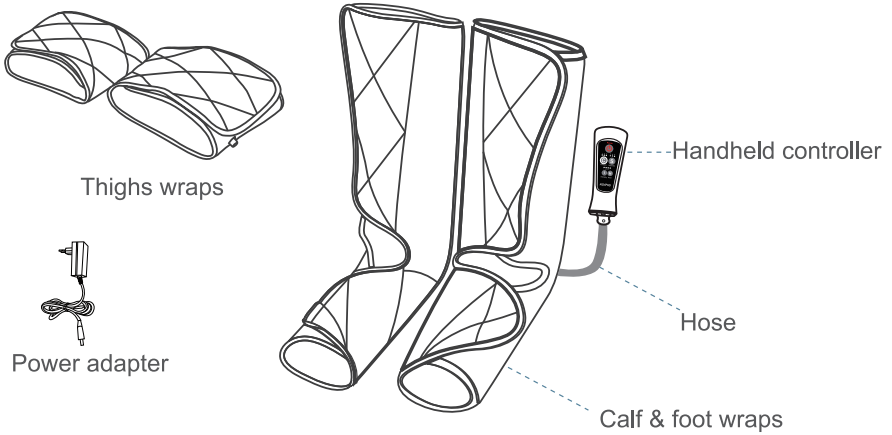
# Using the RENPHO Leg Massager with Heat

## 1. Before Use

1.1 When wearing the wraps, please make sure that you wear the corresponding wrap on the appropriate leg: Left and Right


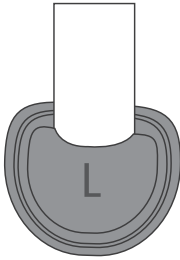
"L" for left leg

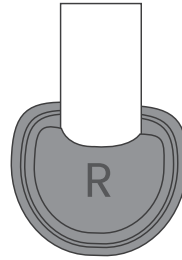"R" for right leg

1.2 Find the mark on the thigh wraps, you can see which side is up.

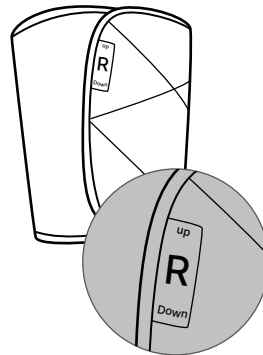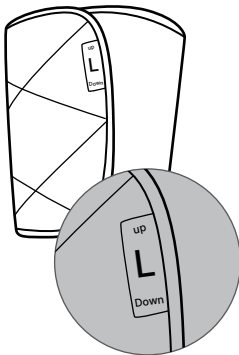
1.3 When you need to use the thigh wrap, make sure the thigh and calf wraps are connected by the air hose. The air hose should be attached securely. Refer to the image below.


1.4 Make sure that the hose has been attached securely by aiming at the snap-fit joint to the controller as shown below.

1.5 If you feel the massage is a bit too powerful or prefer it to be less intense, please follow these recommendations:


- Choose the lowest intensity
- Loosen the velcro-strap of the corresponding part.
- It is recommended to leave enough room for 4 fingers to fit in between the massager and calf. Please refer to the image below.

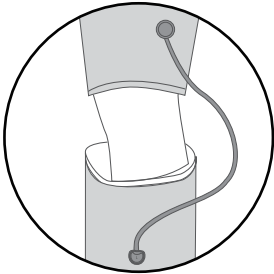1.3

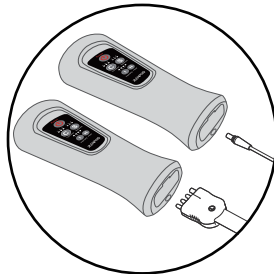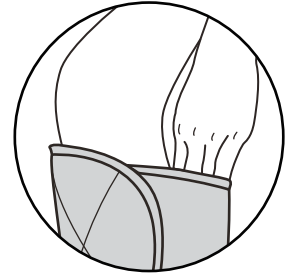
1.4

1.5
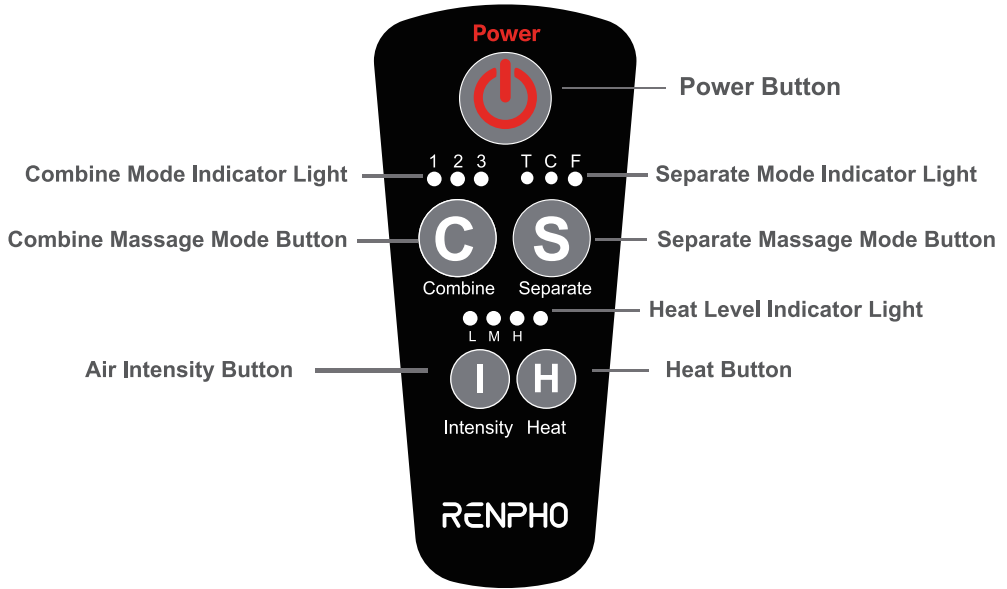
## 2. Operation Instruction

1. “” Power On/Off Button. Press to turn on/off the massager. The default mode is combine message mode 1 with air intensity M. and high heat level(SOLID RED).
2. “” Combine Message Modes. Press to select from the 3 combine message modes.  
Combine 1: Massage full legs.  
Combine 2: Massage feet and calves.  
(It can be used individually, and don't need to connect the air hose of thighs wraps.)  
Combine 3: Massage feet, calves, thighs by turn. It will be turn off when you press the Button C again.

3. “**S**” Separate Massage Modes. Press to select from the thigh, calf and foot separately.  
T: Massages thighs  
C: Massages calves  
F: Massages feet. It will be turn off when you press the Button S again.
4. “**H**” Heat Button: Press to adjust the heat levels (low and high temperature), the default heat level is high temperature. BLINKING RED indicates the low temperature, SOLID RED indicates the high temperature. You can press this button to turn off it.  
Heating area: calves and feet.
5. “**I**” Air Intensity Button. Press to reduce or enhance the intensity. There are 3 intensity levels for each mode, the default intensity is level M. It is recommended to leave a gap available to put 4 fingers in it or choose the lowest intensity if you feel powerful.

**Note:**

1. The massager will be auto shut-off every 20 minutes. And before auto shut off, the indicator light of the corresponding mode will be BLINKING RED. Press the power button again to re-run the massager.
2. The massager comes with memory design, so it automatically operates the last running mode and speed when powering on. (Without power cut off)
3. If you choose one of these functions like combine mode or separate mode or heat, you can turn it off by pressing the button C. Under this situation, when you press the power button, it doesn't come with memory function, the default mode is C1 with air intensity M and high heat level(SOLID RED).

# Maintenance

## 1. Cleaning Method

- **Do not** machine wash. Do not immerse in any liquid to clean. Do not clean it with abrasive cleaners, brushes, gasoline, kerosene, glass polish, or paint thinner. It is recommended to wipe down the massager with a soft, slightly dampened cloth instead.

## 2. Packing Method

Please avoid folding the connection hose when packing; otherwise this may cause inflation failure.

## Troubleshooting

PROBLEM	CAUSE	SOLUTION
The product does not work	The power is not connected or the plug is loose	Plug in the power cord correctly
	The power switch is not on or auto-shut-off timer turns it off	Turn on the power switch
Some part can't inflate	The current massage mode exclude massage this part	Please choose the corresponding mode
	The tube have not been attached securely, and there is air leaks	Please check again if the hose is attached securely to the controller ( no gap between the hose interface and the controller)
	The tube has been folded	<ol style="list-style-type: none"> <li>1. Please avoid folding the tube when packing</li> <li>2. Please contact support@renpho.com</li> </ol>

PROBLEM	CAUSE	SOLUTION
Feel powerful or uncomfortable	Wrap your legs too tight	<ol style="list-style-type: none"> <li>1. Choose the lowest intensity</li> <li>2. Open the Velcro of the corresponding part</li> <li>3. Press POWER button to turn off</li> </ol>
	Maybe use the highest intensity	
The product suddenly stops working	The massager's over-heat protection feature will switch off the massager automatically when too hot	Wait for the massager to cool down before powering it back on
	20 minutes timer turns it off	Press the power button again
The handheld controller gets hot	The product have been used continuously for long time(The temperature rise of handheld control is normal.)	Recommended to use 20 minutes to 30 minutes to avoid OVER USE. After turning off, allow the massager to cool down before using again

PROBLEM	CAUSE	SOLUTION
Feel the difference of intensity when massaging two legs at the same time	When wearing the boots, the gap you leave may be different	Please leave the same gap when you wear it
Can't feel the heat/Heat function is weak	Maybe choose the low temperature	Please choose the high temperature
	Maybe use it just 1-2 minutes	Please use it about 4-6 minutes
If you encounter any other problems, contact: <a href="mailto:support-eu@renpho.com">support-eu@renpho.com</a>		

## Specifications

Product Name: RENPHO Leg Massager with Heat

Model: RP-ALM070H

Product: Input: 12V  $\equiv$  2A

Rated Power: 24 W

AC Adapter: Input: AC 100-240V 50/60Hz

Output: 12V  $\equiv$  2A

Auto-off Timer: 20minutes

## Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

## Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel:** 1800 MY RENPHO (1800 69 7367) (AU ONLY)  
Monday-Friday 9:00AM-6:00PM (AEST)  
+44(0)7434-666088 (UK ONLY)  
Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** [support-au@renpho.com](mailto:support-au@renpho.com) (AU)  
[support-eu@renpho.com](mailto:support-eu@renpho.com)(UK&EU)

\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

### Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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**EU Importer: Renpho GmbH**
**Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.**
**UK REP Richers & Vision Ltd.**
**EU REP C&E Connection E-Commerce (DE) GmbH**