Cordless Massager

Model: RP-EM2020C

CONTENTS

• ENGLISH 01-09

• DEUTSCH 10-19

• FRANÇAIS 20-29

• ITALIANO 30-39

• ESPAÑOL 40-49

Thank you for choosing RENPHO Cordless Massager. Please read the user manual carefully before use and keep it in a safe place for future reference. If you need assistance, please contact our support team.

Caution

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance.

To reduce the risk of electric shock, fire and personal injury, this product must be used as follows:

1. DO NOT reach for an appliance that has fallen into water.
2. DO NOT place in or drop into water or other liquid.
3. DO NOT place or store the device in wet environments or high temperatures.
4. DO NOT drop or improperly use the device, do not remove the screw or attempt to disassemble, do not change or modify the device in any method.
5. DO NOT operate continuously for more than 30 minutes at any given time.
6. Always unplug before cleaning. DO NOT operate or charger the product unattended.

Important Safety Instruction

1. The use or non-rated voltage not (DC 12V \approx 1A) is prohibited.
2. For adults only, it can be used only on dry parts of body and clean skin.
3. Use the massager only in the soft tissues of body without causing pain or discomfort; Do not use on any hard such as bones on the head or the body.
4. Use only massage heads that produce the best results.
5. Whether it is controlled by underground pressure or pressure, it may cause bruises; Often check the massage area, you should stop using immediately once it appears with signs of pain or discomfort.
6. Avoid contact with the massage head and the body at the finger, hair or other parts of body. Beware of pinching.
7. DO NOT press the massage head vigorously during operation.
8. DO NOT use the same part for a long time in a row if you feel unwell, please stop using it immediately. Infants, minors, pregnant women, pacemakers, fractures or metal contents in the body should not use the machine.
9. This product is prohibited under the following conditions:
 - Blood stasis
 - Tumor site
 - Heart disease, with a pacemaker or defibrillatorIf you have experienced joint surgery within 90 days, the product should not be used on the body with 3 inches the surgical position. Sensitive areas: head, face, cervical vertebrae, spine vertebral body. Positions close to the bone with less muscle coverage: the humerus, the back of the foot, the back of the band.
10. DO NOT repair the massager, or replace parts by yourself.
Please contact customer service for any issue.
11. If high temperature or foul odor is experienced while using the


- product, stop using immediately and unplug the device.
- DO NOT put your hands in the gap between the head and the body of the massager when it's working.
 - DO NOT use this product for more than 30 minutes at a time.

About the Renpho Cordless Massager


1. What's in the box

- Cordless Massager x 1
- Massager Head x 5
- Storage Bag x 1
- Adapter x 1
- User Manual x 1

2. Product Introduction

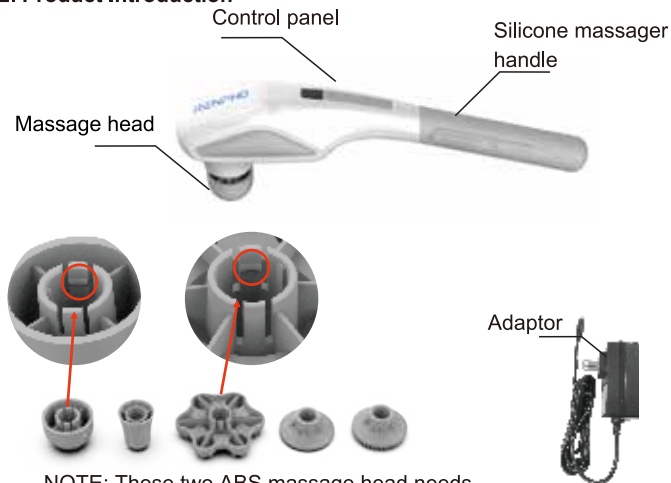NOTE: These two ABS massage head needs to be connected securely by rotating into place at the Snap-Fit Joint.

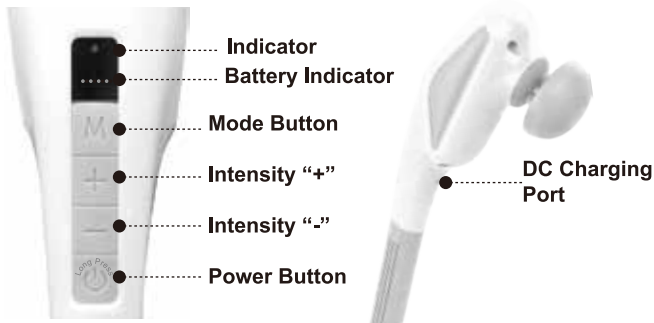Battery Indicator

When the massager is turned on, the indicator will display the corresponding power value through the digital display. (From left to right)

4 green lights on ($\leq 100\%$, $> 75\%$)

3 green lights on ($\leq 75\%$, $> 50\%$)

2 green lights on ($\leq 50\%$, $> 25\%$)

1 green light on ($\leq 25\%$, $> 5\%$)

When the battery power $\leq 10\%$, $> 5\%$, the power indicator will change from a steady green light to a steady red light, until it will stop working.

Head	 Flat massage head	 Intensive massage head	 Round massage head	 Body point massage head	 6-node massage head
Material	TPE	TPE	ABS	TPE	ABS
Applicable body parts	Bulk muscle groups such as arms, thighs	Reflexology & Acupressure & Muscle relaxation	Offer powerful massage	Massage hard-to-reach spots	Deep tissue & Sports massage for waist or back

3. Production Specification

Input: 12V \equiv 1A

Rated Power: 12W

Battery Capacity: 2500mAh

Charging Time: 3 hours

Usage Time: about 140 minutes

Timer: 20 minutes

Product Size: 425x76x127mm/16.7x3.0x5.0in

Net Weight: 1.5kg/3.3lb

Gross Weight: 1.65kg/3.6lb

Operation Instruction

Before use:

Fully charge the massage before first use. The massage runs for about 140 minutes on a full charge. The specific working time depends on how you use it.

1. Plug the adapter into the charging port.
2. During charging, the first green LED indicator on indicates 25% battery power, the the second green LED indicator on indicates 50% battery power, the the third green LED indicator on indicates 75% battery, and the fourth LED indicator stays on green indicates 100% battery power.

Unplug the adapter after full charge, the LED indicator will be off.

3. The battery must be recharged when stored three months.

Noted: The massage is designed to work battery only. It will not work when the AC adapter is connected.

Using the Massager

1. To turn on the massage, long press the power button for 2s. To turn off the massager, short press the power button. The indicator light is always on when using.

The default setting is mode 1 (continuous vibration) and level 3 of intensity.

2. There 5 massage mode settings.

Mode 1: Continuous vibration

Mode 2: Vibration at 1s interval

Mode 3: Vibration at 0.5s interval

Mode 4: Continuous 2s vibration + vibration at 1s interval

Mode 5: Vibration at 0.5s interval*3 + continuous 2s vibration

To change the massage mode, press the Mode button to select a massage mode upon your reference.

3. The massager has 5 intensity settings. To change the intensity, press the "+" to increase the intensity and press the "-" to reduce the intensity.

You will hear a beep with each press. If there is no sound prompt from the buzzer, it means that it is already the highest/lowest level of intensity.

Timing Function: The massager will automatically shut down after 20 minutes of continuous work.

Noted: When you want to change the massage node, please turn off the device firstly.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager doesn't work.	The plug is not well plugged in the socket.	Reinsert the plug.
	Short press the Power button.	Need to long press the Power button to turn on the massager.
Can't hear a tone when pressing the "+"Increase Intensity Button.	The maximum intensity has been reached.	Press the "-" Reduce Intensity Button and you will hear a beep, if not please contact customer support.
Can't hear a tone when pressing the "-"Increase Intensity Button.	The intensity has reached the lowest level.	Press the "+" Increase Intensity Button and you will hear a beep, if not please contact customer support.
The massager shuts off automatically.	The massager will automatically shut down after 20 minutes of continuous work.	Long press the Power button to restart the massager for another 20 minutes.
	The overheat protection function will shut off the massager automatically when it's too hot.	Wait for the massager to cool down before restarting it.
The massage head is sticky.	The heating function is turned on when using 2 additional massage heads.	DO NOT use the heat function when attached the 2 additional massage heads.

Care and Maintenance

1. Place the massager in a safe, cool, and dry place when not in use.
2. Always unplug the massager prior to cleaning.
3. After use, wipe massager with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
4. DO NOT submerge the massager in water or any other liquid.
5. Keep away from all solvents and harsh detergents.
6. DO NOT attempt to repair this massager. There are no user-serviceable parts.
7. DO NOT dispose the massager in a landfill. It contains a Li-ion battery and other electrical components, which are considered hazardous E-waste.

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel:** +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** support-eu@renpho.com(UK&EU)

***For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.**

Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.
EU REP C&E Connection E-Commerce (DE) GmbH
UK REP Richers & Vision Ltd.