RENPHO

User Manual

Foot Massager Compact

RP-FM058R

Black White

CONTENTS

ENGLISH	01 - 07
DEUTSCH	08 - 15
FRANÇAIS	16 - 23
ITALIANO	24 - 31
ESPAÑOL	32 - 39
NEDERLANDS	40 - 47

Thank you for your purchase of RENPHO Foot Massager Compact. The massager offers a powerful massage. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries and/or severe property damage.

Important Safety Instructions

When using this appliance, the following basic precautions should be always followed:

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

1. For household use only, not for medical or commercial purposes.
2. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
3. Children shall not play with the appliance.
4. Cleaning and user maintenance shall not be made by children without supervision.

DANGER – To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. Never use pins or other metallic fasteners with this appliance.
3. Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering or cracking.
4. Do not use or store the appliance in wet environments (e.g. spa, bathroom, sink) or moist conditions, high temperatures or dusty places.
5. Do not reach for the appliance that has fallen into water. Unplug immediately.
6. Do not place or store appliance where it can fall or be pulled into a tub or sink.
7. Do not pour water into the appliance. Do not immerse the appliance in water or other liquids.
8. Do not operate the appliance or handle the plug with wet hands or wet feet.
9. Please attach the plug fully for use to avoid unnecessary electric shock or injury.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

1. If the appliance is to be used by the following user groups, please consult a doctor:
 - 1) Persons suffering from thrombus, severe aneurysm, malignant tumor, venous tumor, skin infections, etc.
 - 2) Persons suffering from osteoporosis, weak bones, abnormal or curved spine, or other acute pains
 - 3) Persons suffering from diabetes or other sensory impairment
 - 4) Persons have a pacemaker or medical devices implanted
 - 5) Pregnant or postpartum women
 - 6) Persons with injured massage areas
 - 7) Persons being treated for a medical condition or physical injury
 - 8) Persons with a heart condition, or abnormal blood circulation on feet
 - 9) Persons with open wounds, discolored areas, or any part of the body that is swollen, burned, inflamed or where skin eruptions or sores are present
 - 10) Persons had received cardiac or visceral surgery
 - 11) Persons insensitive to heat or with heat-feeling loss
 - 12) Persons suffering from a fever
 - 13) Persons feeling discomfort
 - 14) Persons with carsick, seasick or reacting vigorously to other means of transport
 - 15) Persons easy to get sprained or injured

Otherwise it may cause accidents, injuries or discomfort!

If you have any discomfort (such as inflamed skin, redness, itching, etc.), please stop using the appliance immediately and consult your doctor.

2. The appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before cleaning, putting on or taking off parts.
3. Do not use the appliance while covered in quilt or blanket or near heat resources (e.g. radiators or stoves). Excessive heating can occur and cause fire, electric shock, or injury to persons.
4. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
5. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
6. Do not carry this appliance by supply cord or use cord as a handle. Keep the cord away from heated surfaces.
7. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
8. Never drop or insert any object into any opening.
9. Do not use outdoors. Place and use the appliance on a dry, hard and flat surface indoors.
10. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
11. Do not crush – Avoid sharp folds. Keep away from sharp objects to avoid scratching.
12. Do not stand on or in appliance. Use only while seated.
13. Do not place excessively heavy objects on the appliance to avoid damage.
14. Do not allow pets and animals to lie or chew on this appliance or on the electric cord.
15. Do not use while sleeping/feeling asleep or after intake of alcohol.
16. Any modification or disassembly of the appliance without authorization is prohibited.
17. If a high temperature or foul odor occurs, please stop using the appliance immediately and cut off power.
18. Do not put your feet reversely into the appliance, otherwise the malfunction may occur.
19. Make sure the massage is completed before removing your feet.
20. To disconnect, turn all controls to the off position, then remove plug from outlet.
21. Only use the original adapter provided by RENPHO.
22. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
23. Use heated surfaces carefully. May cause serious burns. Do not use it on insensitive skin or on a person with poor blood circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
24. Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
25. The appliance is only to be used with the power supply unit provided with the appliance.
26. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
27. The appliance must not be located immediately below a socket-outlet.
28. Do not operate this appliance if the cord or plug is damaged, if it is not functioning properly, if it has been dropped or damaged, or if it has been dropped into water. Please return the appliance to a service center for examination and repair.
29. Do not use the machine while under the influence of alcohol, sleeping pills or sedatives.

About the Battery

1. Do not expose the appliance or battery to excessive temperatures.
2. Exhausted batteries are to be removed from the appliance and safely disposed of.
3. If the appliance is to be stored unused for a long period, the batteries should be removed.
4. Do not use modified or damaged batteries.
5. Be aware of the risk of terminals of the battery-operated appliance or

battery being short-circuited by metal objects.

6. This appliance contains non-rechargeable batteries, these batteries are not to be recharged.
7. Different types of batteries or new and used batteries are not to be mixed.
8. Remove and immediately recycle or dispose of used batteries according to local regulations and keep away from children.
9. Do NOT dispose of batteries in household trash or incinerate.
10. Even used batteries may cause severe injury or death.
11. Call a local poison control center for treatment information.
12. Clean the battery contacts and also those of the equipment prior to battery installation.
13. Do not force discharge, recharge, disassemble, heat above (manufacturer's specified temperature rating) or incinerate. Doing so may result in injury due to venting, leakage or explosion resulting in chemical burns.
14. Ensure the batteries are installed correctly according to polarity (+ and -).
15. Do not mix old and new batteries, different brands or types of batteries, such as alkaline, carbon-zinc, or rechargeable batteries.
16. Remove and immediately recycle or dispose of batteries from equipment not used for an extended period of time according to local regulations.
17. Always completely secure the battery compartment. If the battery compartment does not close securely, stop using the product, remove the batteries, and keep them away from children.
18. Always select the correct size and type of battery most suitable for the intended use. Information provided with the equipment to assist correct battery selection should be retained for reference.
19. In case of ingestion of a cell or battery, seek medical assistance promptly. Swallowing lithium coin cells or batteries can cause chemical burns, perforation of soft tissue, and in severe cases can cause death. They must be removed immediately if swallowed.
20. Do not allow children to replace batteries without adult supervision.
21. Do not force discharge batteries. Do not deform batteries. Do not dispose of batteries in fire.
22. Exhausted batteries should be immediately removed from equipment and properly disposed of. When discharged batteries are kept in the equipment for a long time, electrolyte leakage can occur causing damage to the equipment and/or personal injury.
23. A lithium battery with a damaged container should not be exposed to water.
24. Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.

About the RENPHO Foot Massager Compact

1. What's in the box

- 1 x Foot Massager Compact
- 1 x Remote Control
- 1 x Adapter
- 1 x User Manual

2. Specifications

- Product Name: Foot Massager
- Model: RP-FM058R
- Input: 12V 2.5A
- Rated Power: 30W
- Auto-off Timer: 15/20/25/30min
- Heating Area: Bottom
- Material: ABS + Lycra Fabric
- Heating Temperature: Below 50 °C/122 °F
- Product Weight: 3.4kg/7.5lb

3. Product Overview

Detachable foot cover with zipper

Touch panel

Remote control

Power


Time (Auto-off)

Air intensity

Knead intensity

Heat

4. Function

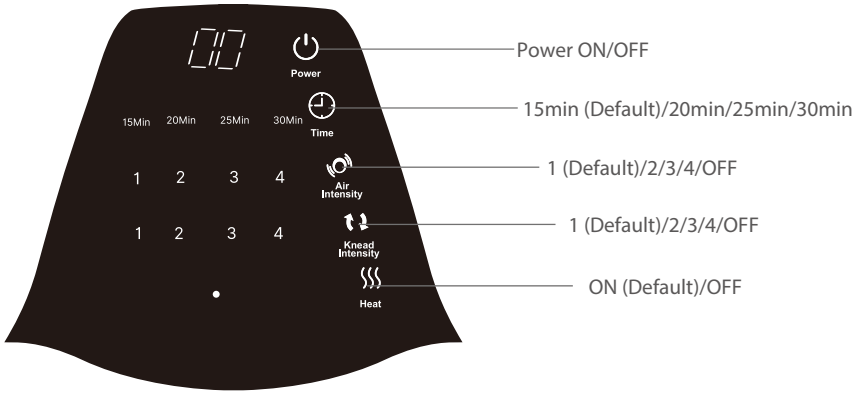Using the RENPHO Foot Massager Compact

It is recommended to enjoy the massage for no more than 30 minutes at a time. Longer massage time may cause mild pain.

Connect the plug to the power source. The Power Button flashes green accompanied by a beep sound and the massager enters standby mode. Remove the protective film from the touch panel before use.

1. Press the Power Button to turn on the massager (Power Button will stay on green while other buttons on white). Check the default time and intensity in the above Function part. Press again to enter the standby mode (Only Power Button will flash green).
2. Press the Time Button to switch the auto-off time among 15, 20, 25, or 30 minutes. The timer will reset each time you press the Time Button.
3. Press the Air Intensity Button to adjust the air intensity among 1/2/3/4/OFF.
4. Press the Knead Intensity Button to adjust the kneading intensity among 1/2/3/4/OFF.
5. Press the Heat Button to turn on/off the heating function.

Note:

1. The air pressure and kneading function can work independently. The heating function can work only when the air pressure or kneading function is on.
2. Each press is accompanied by a beep sound and corresponding indicator light (When off, the LED goes out).
3. If none of the heating, air pressure and kneading function is working, and the Time and Power button indicator lights are on, the massager will enter the standby mode if there is no further operation within 10 seconds.
4. The massager is programmed with a memory function. It will operate in the pattern set last time once you turn on the Power Button again. The memory function will not work if only the Timer function is working before turning off the massager, or if you cut off the power supply.
5. The massager fits up to US Men's Size 12 (EU Men's size 46).
6. Detachable foot cover: unzip to remove and clean.

Product Care and Maintenance

1. Unplug the appliance and allow it to cool before cleaning. You can use a soft cotton cloth dampened with neutral detergent to wipe it clean.
2. The use of non-rated voltage not (DC 12V, 2.5A) is prohibited.
3. It is prohibited for non-professional personnel to repair the appliance or replace parts.
4. If the power adapter is damaged, it must be replaced by the maintenance department or qualified personnel approved by the manufacturer so as to avoid potential hazards.
5. Do not use or store the appliance in wet environments or high temperatures.
6. Do not use sharp objects to handle the appliance in order to avoid damage.

Troubleshooting Guide

Problem	Cause	Solution
The massager does not work	<ol style="list-style-type: none">1. Disconnected power or loose plug causing poor contact.2. The power button is off or the timer has been reset.3. The motor auto shuts-off due to the device being on too long activating overheating protection.	<ol style="list-style-type: none">1. Re-plug to the power source.2. Turn on the power button.3. Use the massager after the motor cools off naturally.
The massage intensity is too high or too low	<ol style="list-style-type: none">1. Your foot is positioned at the wrong angle.2. Didn't choose intensity level which is right for you.	<ol style="list-style-type: none">1. Adjust the contact angle to get moderate intensity.2. The massager has 4 intensity levels, choose the intensity you prefer. We recommend low intensity to begin with. Or individually choose air pressure or kneading function to reduce the intensity.
No heating or insufficient heating	<ol style="list-style-type: none">1. The heating function is turned off.2. The product structure and working principle limits how fast the bottom can heat up. Therefore, the heating-up may seem a bit slow, 7-10 minutes before you can feel the obvious heat.	<ol style="list-style-type: none">1. Turn on the Heat Button.2. It takes about 7-10 minutes to fully heat up. Please be patient as this is to avoid the temperature rising too fast.
The touch panel is not working properly	<ol style="list-style-type: none">1. The plastic film on the touch panel is not removed.2. Fingers are wet.3. Insufficient touch area.	<ol style="list-style-type: none">1. Remove the plastic film before use.2. Dry your fingers before touching the button.3. Use your thumb to press buttons to get better touch sensitivity.

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:
<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

-  Tel: +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)
-  Email: support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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