
# RENPHO

Smart Health Simplified

## User Manual

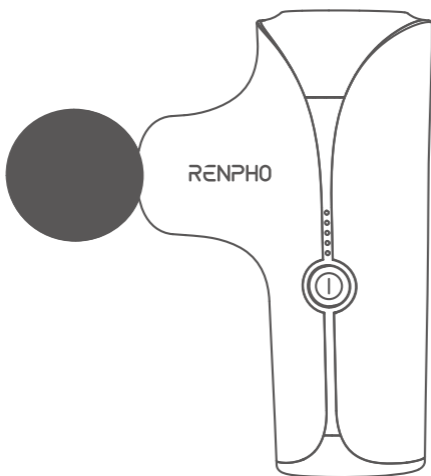Model: RP-GM173

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# CONTENTS

- 
- ENGLISH 01-15

- 
- DEUTSCH 16-30

- 
- FRANÇAIS 31-44

- 
- ITALIANO 45-58

- 
- ESPAÑOL 59-72
-

Thank you for purchasing RENPHO mini massage gun. The massage gun offers a powerful percussion massage. While enjoying muscle and body relief, you should use it with care. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

## **Introduction**

Mini Massage Gun is a device for deep muscle and fascia relaxation. It is suitable for physical massage after exercise. It can reduce muscle stiffness and pain, alleviate muscle fatigue and accelerate blood circulation to achieve the purpose of improving the health of body soft tissue.

## **Important Safety Instructions**

### **READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.**

When using this device, the following basic precautions should be always followed:

1. The device is intended for household use only and not intended for medical use.
2. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the device.

### **DANGER - To reduce the risk of electric shock:**

1. Do not attempt to retrieve the device if it falls into water.
2. Do not use while bathing or in a shower.
3. Do not use the device with wet hands.

4. Do not place or store the device in areas where it may fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

**WARNING - To reduce the risk of fire, electric shock, or injury to persons:**

**1. Please consult your doctor before using if you are a user of the following persons:**

- 1) Children under 8
- 2) Pregnant women
- 3) Persons with blood stasis
- 4) Persons with fractures
- 5) Persons with heart disease, a pacemaker, defibrillator, or metal contents in the body
- 6) If you have experienced joint surgery within 90 days, do not use the device on your body within 3 inches of the surgical site.

**If you feel unwell or have any unusual discomfort, please stop using the device immediately.**

**2. Do not use the massage gun on**

- 1) Sensitive skin, swollen or inflamed areas, in areas where there is poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain
- 2) Head, face, cervical vertebrae, body bones, varicose veins, tumor site, throat area, any possible fracture site, positions close to the bony areas: the humerus, the back of the foot, the back of the hand, etc
- 3) Genital area (Do not use this device for sexual purposes. It is a powerful massaging device and may cause severe bodily harm.)
- 4) Unconscious or sleeping persons
3. Do not operate under a blanket or pillow. Excessive heating may occur and cause fire, electric shock, or injury to persons.
4. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
5. Close supervision is necessary when this device is used by, on or near children, invalids, and disabled persons.

6. Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
7. Never operate this device if it is not working properly, has been damaged, or dropped into water.
8. Never drop or insert any object into any opening.
9. Do not operate in the areas where aerosol (spray) products are being used or where oxygen is being administered.
10. Do not keep the device near a magnetic field or in the areas where it can be affected by a magnetic field as it may cause batteries to discharge.
11. Avoid touching the gap between the massage head and body with your finger or hair to avoid pinching.
12. Keep the device away from heat and avoid using it in places where the surrounding temperature is relatively hot.
13. Do not drop, apply excessive force or heavy weight over the device to avoid malfunction and/or damage.
14. If a high temperature or foul odor occurs while using the device, stop using the device immediately.
15. Do not use the device in close proximity to long hair, loose clothing or jewelry.
16. Avoid applying excessive pressure on the massage head and refrain from massaging the same area continuously for an extended period of time.

## **About the Battery**

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.

6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. **CAUTION** - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

## SAVE THESE INSTRUCTIONS

### About the RENPHO Mini Massage Gun

#### 1. What's in the box

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Massage Head x 4

Storage Case x 1


Type-C Cable x 1

Strap x 1

User Manual x 1

## 2. Product Introduction

Removable and replaceable massage head

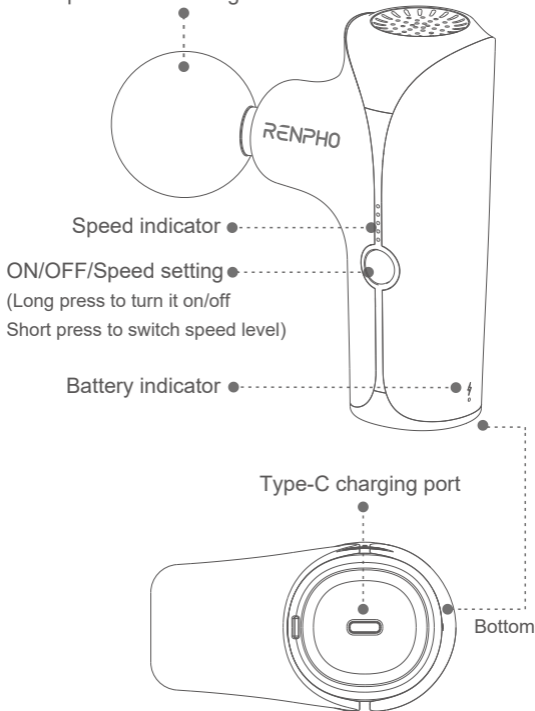EN-05

Damping  
Head

U-shaped  
Head

Ball Head

Bullet Head

### 3. Specification

Product Name: Mini Massage Gun

Model: RP-GM173

Input: 5V  $\overline{\text{---}}$  2A

Battery Voltage: DC 11.1V

Battery Capacity: 2500mAh/27.75Wh

Amplitude: 8mm

Vertical Vibration Frequency  $\leq$  60Hz

Speed Setting: 1-5 levels 1800-3200r/min

Charging Time: 3 hours

Operating Time: About 4 hours

Auto-off Timer: 10 minutes

Product Dimensions: 139 x 87 x 49mm

5.5 x 3.4 x 1.9inch

Product Weight: 0.45kg/0.99lb

Gross Weight: 0.9kg/1.98lb

# Using Mini Massage Gun

## 1. Important tips before use

### 1.1 Charging

● Please fully charge the massage gun for about 3 hours before use.

● Connect the output end of the connector to the DC charging port at the bottom of the massager, and then connect the input end of the connector to the mobile power supply or the mobile phone charger. (5V = 2A, if you use less than 2A charger, the time will be longer.)

● When the battery power is more than 30%, the battery indicator stays Green; When the battery power is less than 20%, the indicator light is flashing in red; When the power is less than 5%, please charged it immediately, otherwise it will automatically shutdown after flashing in red for few seconds. It is not recommended that the battery power be completely depleted to 10% level.

● For safety considerations, the massage gun will NOT work while charging. Additionally, the massage gun will automatically shut off after 10mins of continuous use. If you want longer massage, please turn it on then.

### 1.2 Battery Life

● This Massage gun has a built-in rechargeable battery. A full charge takes about 3 hours. Total battery life is up to about 4 hours after a full charge.


**DON'T** use the massage gun when charging!

● Battery life can vary by application and/or operating conditions. If the product is not in use for a long time, please keep it powered off and charge it once a month so that the battery will last longer.

● Please keep the massage gun away from heat and avoid using it in places where the surrounding temperature is relatively hot.

## 2. Operation Instruction

2.1 Select a preferred massage head and screw it fully into the slot of the device until it is securely attached. Ensure the device is turned off when installing the massage head.

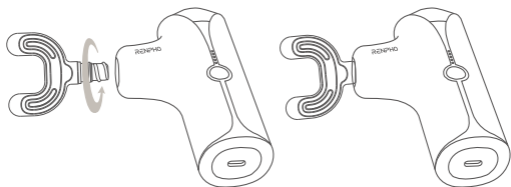2.2 Turn on the device by long pressing the power button. Short press the power button, the gear will be switched between 1-5 levels. If the power button is long pressed for 2 seconds, the device will stop immediately. LED indicator indicates the corresponding level, the lowest level of 1, the highest level of 5.

2.3 Apply the required pressure with the preset vibration massage frequency to massage the corresponding body parts.


\* After 10mins of operation, the device automatically will power off.

Shoulder Massage


Back Massage

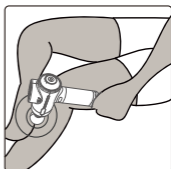Foot Massage

Arm Massage

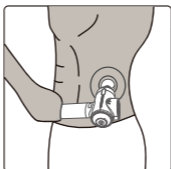Waist Massage

Leg Massage

## Product Care and Maintenance

- Place the massage gun in a safe, cool, and dry place when not in use.
- Always unplug the massage gun prior to cleaning.
- After use, wipe massage gun with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
- DO NOT submerge the massage gun in water or any other liquid.
- Keep away from all solvents and harsh detergents.
- DO NOT attempt to repair this massage gun. There are no user- service-able parts.
- DO NOT dispose the massage gun in a land fill. It contains a Li-ion battery and other electrical components, which is considered hazardous E-waste.
- Do not disassemble the screw or attempt to disassemble the equipment.
- When the equipment is not used for a long time, the battery should be charged once a month.
- Massage equipment should be stored in a well ventilated, non-corrosive gas warehouse, 250 mm off the ground, no more than 5 layers of massager stacked.

## Environmental conditions

Because the electromagnetic field produced by massage motor is relatively small, it will not affect the electronic equipment and human body in general! If the main engine (including lithium batteries) and accessories of this instrument reach the service life, please do not discard them at will. Otherwise, electronic waste will pollute the environment. Please make corresponding treatment according to the relevant regulations of local environmental protection department for the treatment of such electronic waste.

Temperature:  $-40^{\circ}\text{C} \sim +50^{\circ}\text{C}$  /  $-104 - 122^{\circ}\text{F}$

Relative humidity:  $\leq 90\%$

Barometric pressure:  $70\text{kPa} \sim 106\text{kPa}$

## Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	Power is off.	Long press the Power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically.	The massage gun's overheat protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before powering it back on.
	Massage gun shuts off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is too inserted too tight.	Plug out the massage head a little, then power it on.

## Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:  
<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

# Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **TEL: +44(0)7434 666088 (UK ONLY)**  
Monday-Friday 9:00AM-5:00PM (GMT)  
**+1(844) 417 0149 (US&CA)**  
Monday-Friday 9:00AM-4:30PM  
**1800 MY RENPHO (1800 69 7367) (AU ONLY)**  
Monday-Friday 9:00AM-6:00PM (AEST)  
**+49 40 23969285 (DE ONLY)**  
Monday-Friday 9:00AM-6:00PM (CET)

 **Email: support-eu@renpho.com (UK&EU)**  
**support@renpho.com (US&CA)**  
**support-au@renpho.com (AU)**

\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

# FCC Regulatory Compliance

**Warning:** Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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### 1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

### 2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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