RENPHO

User Manual

Shoulder & Neck Massager

Model: RP-SNM061

CONTENTS

1. ENGLISH	01-07
2. DEUTSCH	08-15
3. FRANÇAIS	16-22
4. ITALIANO	23-29
5. ESPAÑOL	30-36

Thank you for your purchase of RENPHO Kneading Shoulder Massager. The massager offers a powerful shiatsu massage. While enjoying massage relief, you should use it with care and safety in mind. To ensure safety and avoid any property damage, please read this manual carefully and use it as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

Important Safety Instructions

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the device.
3. This device is not intended for medical use. It just offers a massage function.
4. Please consult a doctor before use if you fit into any of the following user groups:
 - Persons who are pregnant or postpartum
 - Persons with skin infections, allergies, redness and/or skin burns
 - Persons with osteoporosis
 - Persons who had undergone surgery
 - Persons with any swelling or wound around the abdomen
 - Persons with a pacemaker or other implantable medical devices
 - Patients with heart disease, malignant tumors, or other serious health conditions
 - Persons who are receiving medical treatment
 - Persons with heart problems, abnormal blood pressure, diabetes or neurological disorders
5. If you feel unwell or in any unusual discomfort, please stop using the appliance immediately.
6. NOT suitable for children.
7. DO NOT use it while driving.
8. The elderly need to use under the doctor's guidance.
9. DO NOT use this appliance after drinking or taking sleeping pills.
10. NEVER try to inspect or open the appliance in any way.
11. DO NOT use this appliance for too long, as this may cause minimal damage to the skin.
12. DO NOT put heavy objects on the appliance.
13. NEVER use pins or other metallic fasteners with this appliance.
14. DO NOT use this appliance if you apply any oils or lotions to your skin.
15. NEVER clean this appliance unless it is unplugged. Use a cloth dampened with a very small amount of water to clean. You can use non-abrasive cleaners. DO NOT plug this appliance back in until it is totally dry.

CAUTION - To reduce the risk of burns:

1. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.

2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
4. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
5. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
6. Do not carry this appliance by supply cord or use cord as a handle.
7. Keep the cord away from heated surfaces.
8. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and etc.
9. Never drop or insert any object into any opening.
10. Do not use outdoors.
11. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
12. To disconnect, turn all controls to the off position, then remove plug from outlet.
13. The appliance is only to be used with the power supply unit provided with the appliance.
14. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.


SAVE THESE INSTRUCTIONS

About the RENPHO Shoulder & Neck Massager

1. What's in the box

- Shoulder & Neck Massager x 1
- User manual x 1
- Power adapter x 1
- Car adapter x 1
- Soft sleeve x 1


2. Product overview

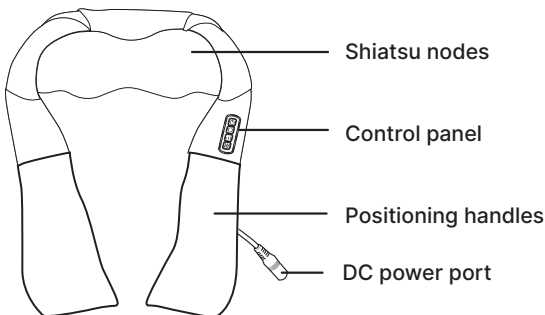Power adapter


Soft sleeve

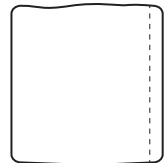Car adapter

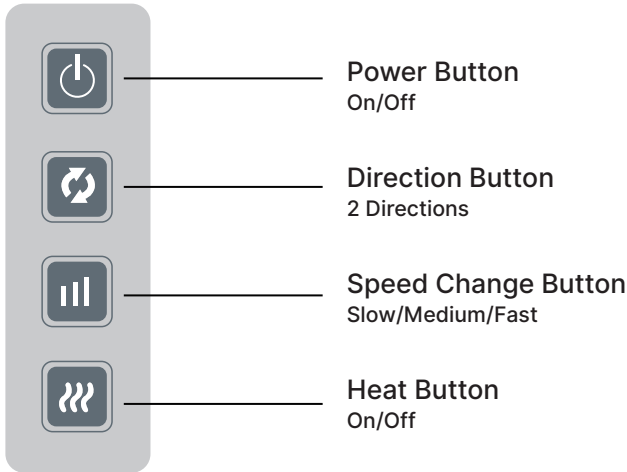Using the RENPHO Shoulder & Neck Massager


1. Connect the power adapter into the DC power port, then plug the power adapter into an electric outlet. (Car adapter included. You can also use it in your car.)
2. Grip the positioning handles and drape the main unit around the area you want to massage.
3. Press  the Power button. The massager will turn on. The default mode is as below:
Medium speed
Heat on
15-minute timer default
(The direction of rotation changes every minute)
4. Press  the Direction button to change the massage direction.
5. Press  the Speed Change button to change the speed.
6. Press  the Heat button to turn off or on the heat function.
7. Use the positioning handles to reposition the Shiatsu Nodes to the area you want to massage.
8. Press  the Power button at any time to turn the appliance off.
9. The appliance shuts off automatically after 15 minutes of continuous use.


Massage Areas


Neck&Shoulder

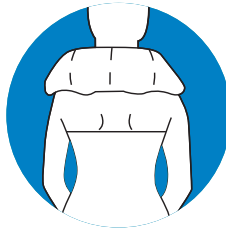
Neck

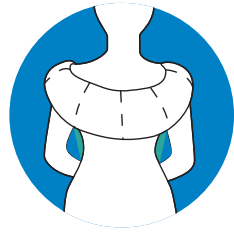Back

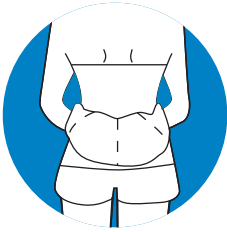Waist

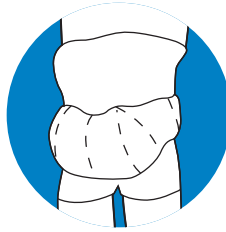Abdomen


Thighs

Legs

Feet

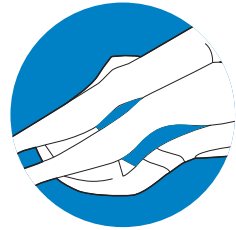Calves

Product Care and Maintenance

1. Place the appliance in a safe, cool, and dry place when not in use.
2. After unplugging, wipe with a soft, damp cloth using non-abrasive cleaners.
3. Never submerge the appliance in any liquid.
4. Keep away from all solvents and harsh detergents.
5. **DO NOT** attempt to repair this appliance. There are no user-serviceable parts.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager doesn't turn on	Power adapter plug is loose	Reconnect the plug to the outlet
	Power is off	Press the power button to turn on the massager
	No power	Try a different power outlet
Shuts off automatically	The 15-minute timer	The massager has a timer and will turn off on its own after 15 minutes. You can press the power button to start another 15-minute massage
	The massager's overheat protection feature will switch off the massager automatically when too hot	Wait for the massager to cool down before powering it back on
Can't feel heat	You used the soft sleeve	Take off soft sleeve
	To prevent overheating, our temperature is controlled at 118°F	It's easy to hurt your neck when the heat is too strong
Feel pain while using	The intensity of the deep kneading massage is strong	Decrease speed
		Add the soft sleeve on the massager to further reduce the intensity of the deep kneading massage
Can't feel intensity	Your pull on the handle is not strong enough	Increase your pull on the handle
	The intensity of the deep kneading massage is weak	Increase speed

Technical Specifications

Product Name	Shoulder & Neck Massager
Model	RP-SNM061
Rated Voltage	DC 12V
Rated Power	24W
Standby Power	<0.5W
Noise	≤55dB

Technical Specifications

Temperature	< 118°F
Colour	Gray
Gross Weight	6.17lb
Size	18.89(L) × 6.70(W) × 7.87(H) inch
Default Timer	15 Minutes
What's Inside	Shoulder & Neck Massager, Power Adapter, Car Adapter, User Manual, Soft Sleeve

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not void the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a retail store or other collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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