# RENPHO

Model: RP-SNM061

## Kneading Shoulder Massager

# CONTENTS

---

1. ENGLISH .....	1-6
2. DEUTSCH .....	7-13
3. FRANÇAIS .....	14-20
4. ITALIANO .....	21-27
5. ESPAÑOL .....	28-34

Thank you for your purchase of RENPHO Kneading Shoulder Massager. The massager offers a powerful shiatsu massage. While enjoying massage relief, you should use it with care and safety in mind. To ensure safety and avoid any property damage, please read this manual carefully and use it as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

## Important Safety Instructions

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance.

Cleaning and user maintenance shall not be made by children without supervision.



- This product is not intended for medical use. It just offers a massage function.
- Any Individual who may be pregnant or recently gave birth, has body diseases or skin diseases, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor's care should consult with a medical professional before using a massaging device designed for home use.
- If you feel unwell or in any unusual discomfort, please stop using the product immediately.
- **NOT** suitable for children.
- **DO NOT** use it while driving.
- The elderly need to use under the doctor's guidance.
- **DO NOT** use this product after drinking or taking sleeping pills.
- **NEVER** try to inspect or open the product in any way.
- **DO NOT** use this product for too long (Which may cause minimal damage to the skin).
- **DO NOT** put heavy objects on the product.
- **DO NOT** continue to use this product after it has been damaged.
- **DO NOT** drop the product on the ground at high altitude.
- **DO NOT** reach for an appliance that has fallen into the water. Unplug immediately.
- **DO NOT** place in or drop into water or other liquids.
- **NEVER** use pins or other metallic fasteners with this appliance.
- **DO NOT** use this product if you apply any oils or lotions to your skin.
- **DO NOT** carry this product by the power cord.
- **DO NOT** allow the power cord near any heated surfaces.
- **NEVER** clean this product unless it is unplugged. Use a cloth dampened with a very small amount of water to clean. You can use non-abrasive cleaners. **DO NOT** plug this product back in until it is totally dry.
- **DO NOT** operate under a blanket or pillow. Excessive heating can occur or cause fire, electric shock or injury to persons.
- Power off and disconnect products from the mains after use.

# About the RENPHO Neck and Shoulder Massager

## 1. What's in the box

- Kneading shoulder massager x 1
- User manual x 1
- Power adapter x 1
- Car adapter x 1
- Soft sleeve x 1


## 2. Production Instruction

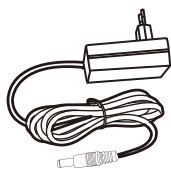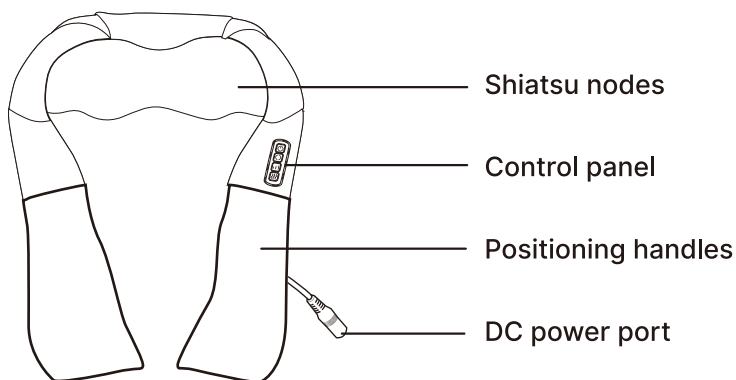Power adapter


Car adapter

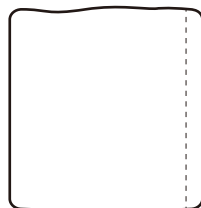Soft sleeve
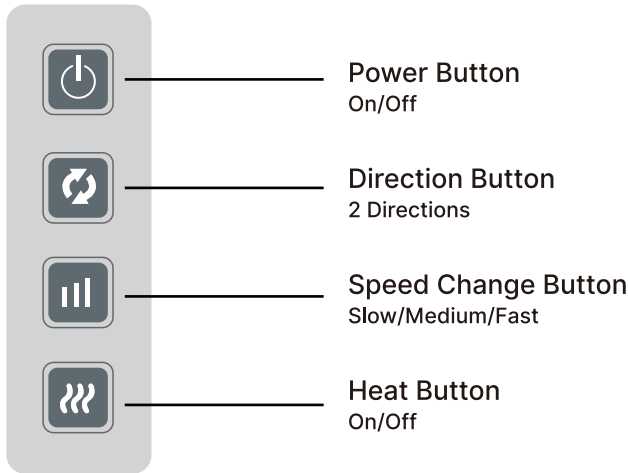
## Using the Neck and Shoulder Massager


1. Connect the power adapter into the DC power port, then plug the power adapter into an electric outlet. (Car adapter included. You can also use in your car)
2. Grip the positioning handles and drape the main unit around the area you want to massage.
3. Press  the Power button. The massager will turn on. The default mode is as below:  
Medium speed  
Heat on  
15 min default  
(The direction of rotation changes every minute)
4. Press  the Direction button to change the massage direction.
5. Press  the Speed change button to change the speed.
6. Press  the Heat button to power off or on heat function.
7. Use the positioning handles to reposition the Shiatsu Nodes to the area you want to massage.
8. Press  the Power button at any time to turn the unit off.
9. The unit shuts off automatically after 15 minutes of continuous use.


## Massage Areas


Neck&Shoulder

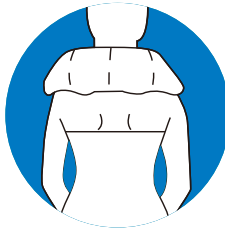
Neck

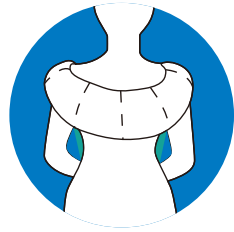Back

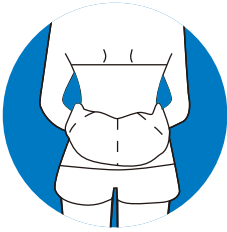Waist

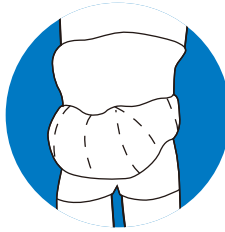Abdomen


Thighs

Legs

Foot

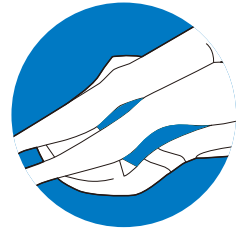Calf

## Product Care and Maintenance

1. Place the massager in a safe, cool, and dry place when not in use.
2. After unplugging, wipe with a soft, damp cloth using non-abrasive cleaners.
3. Never submerge the unit in any liquid.
4. Keep away from all solvents and harsh detergents.
5. **DO NOT** attempt to repair this massager. There are no user-serviceable parts.

## Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager doesn't turn on	Power adapter plug is loose	Reconnect the plug to the outlet
	Power is off	Press the power button to turn on the massager
	No power	Try a different power outlet
Shuts off automatically	The 15-minute timer	The massager has a timer and will turn off on its own after 15 minutes. You can press the power button to start another 15 minute massage
	The massager's overheat protection feature will switch off the massager automatically when too hot	Wait for the massager to cool down before powering it back on
Can't feel heat	You used the soft sleeve	Take off soft sleeve
	To prevent overheating, our temperature is controlled at 118°F	It's easy to hurt your neck when the heat is too strong
Feel pain while using	The intensity of the deep kneading massage is strong	Decrease speed
		Add the soft sleeve on the massager to further reduce the intensity of the deep kneading massage
Can't feel intensity	Your pull on the handle is not strong enough	Increase your pull on the handle
	The intensity of the deep kneading massage is weak	Increase speed

## Specifications

Product Name	RENPHO Kneading Shoulder Massager
Model	RP-SNM061
Rated Voltage	DC 12V
Rated Power	24W
Standby Power	< 0.5W
Noise	≤55dB

## Specifications

Temperature	< 118°F
Colour	Grey
Gross Weight	6.17lb
Size	18.89(L) × 6.70(W) × 7.87(H) inch
Default Timer	15 Minutes
What's Inside	Massager, Power Adapter, Car Adpater, User Manual,Soft Sleeve

## Warranty (1-Year Period)

The manufacturer's warranty is void under the following circumstances but not limited to :

- Improper or inadequate maintenance or modification.
- Accident, misuse, abuse, contamination, or other external causes.
- The use of unspecified adapter and accessories.
- Loss or damage in transit.
- Damages that occur as a result of one's failure to follow the instructions.
- This warranty does not apply to expendable or consumable parts and does not extend to any product from which the serial number has been removed.

## Customer Service

1. If you have any product questions or concerns, don't hesitate to contact us directly at [support@renpho.com](mailto:support@renpho.com) or call our customer service hotline:

-  **Tel:** +44(0)7434-666088 (UK ONLY)  
Monday-Friday 9:00AM-5:00PM(GMT)
-  **Email:** [support-eu@renpho.com](mailto:support-eu@renpho.com)(UK&EU)
-  **Web:** [www.renpho.com](http://www.renpho.com)

2. We have an award-winning warranty, exchange, and customer service team that guarantees hassle-free solutions to any issue you might have within the 24 business hours.

### Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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## Declaration of Conformity

We (EU Representative is responsible for this declaration)

C&E Connection E-Commerce (DE) GmbH  
(Company name)

Zum Linnegraben 20, 65933, Frankfurt am Main, Germany  
(Company address)

declare under our sole responsibility that the product

Tradename: RENPHO

Equipment: Kneading Shoulder Massager

Model No. : RP-SNM061

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to EMC Directive(2014/30/EU) &RoHS(2002/95/EC) product is responsible to affix **CE** marking, the following standards were applied:

- EN 55014-1:2006+A1:2009+A2:2011*
- EN 55014-2:2015*
- EN 61000-3-2:2014*
- EN 61000-3-3:2013*
- IEC 62321 Series*

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Datum: May.20.2020

# Declaration of Conformity

We (UK Representative is responsible for this declaration)

Richers & Vision Ltd.  
(Company name)

Unit 58B, Bordon Enterprise Park, Budds Lane, Bordon, Hampshire, GU35 0JE, U.K.  
(Company address)

declare under our sole responsibility that the product

Tradename: RENPHO

Equipment: Kneading Shoulder Massager

Model No. : RP-SNM061

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to EMC Directive(2014/30/EU) &RoHS(2002/95/EC) product is responsible to affix **CE** marking, the following standards were applied:

*EN 55014-1:2006+A1:2009+A2:2011*

*EN 55014-2:2015*

*EN 61000-3-2:2014*

*EN 61000-3-3:2013*

*IEC 62321 Series*

Full Name: Comy Xiao

Position: Approbation Specialist

Signature: 

Date: May.20.2020

**EU|REP C&E Connection E-Commerce (DE) GmbH**

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Unit 58B, Bordon Enterprise Park, Budds Lane, Bordon, Hampshire, GU35 0JE, U.K.

**Shenzhen Runfeng Digital Technology Co., LTD**

Chuanghui Building, Intersection of Banxuegang Boulevard and Yongxiang Road, Bantian Subdistrict, Longgang District, Shenzhen

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**UK|REP Richers & Vision Ltd.**